

PEEK-OPTIMA™ POLYMER KNEE – A More Natural Solution to Total Knee Replacement

AUTHOR: Ian Revie, PhD – Invibio Biomaterial Solutions

Total Knee Replacement Satisfaction Levels

According to the Organization for Economic and Co-operation Development (OECD) health statistics, in 2015, there were approximately 2.6 million total knee replacement (TKR) surgeries reported for the 35 member countries.¹ Yet, a significant number of people who have undergone TKR are unhappy with the results.² An estimated 20% of TKR patients, totaling about 520,000 worldwide, are unsatisfied.²⁻³


Why are Patients Dissatisfied?

The literature proposes various reasons for TKR dissatisfaction, but little is proven.⁴⁻⁵ One researcher linked positive patient expectation to improved outcomes.⁴ Moreover, fulfilling kneeling, squatting and stair climbing expectations correlated highly with satisfaction.⁴ In fact, preparing patients for surgery and setting an expectation of recovery timeframes and outcomes resulted in a 5% improvement in patient satisfaction.⁴

Surgical precision and materials used may also play a part in patient satisfaction levels. Some unsatisfactory TKR results were due to primary malalignment of the femoral component.⁵ Different femoral component materials can also affect a patient's ability to perform normal daily activities post-operatively. TKR patient feedback across social media and forums listed TKR post-operative side effects ranging from knees feeling heavy or loose, to being cold in the winter, to knee creaking, popping or rattling, which may all potentially be associated with metal.

Introducing the PEEK-OPTIMA Polymer Knee Implant

The PEEK-OPTIMA Knee Implant has been developed with these problems in mind. It offers surgeons and patients a femoral component that has the potential to improve both operative procedure and post-operative outcome. Most importantly the goal is to improve a patient's quality of life.

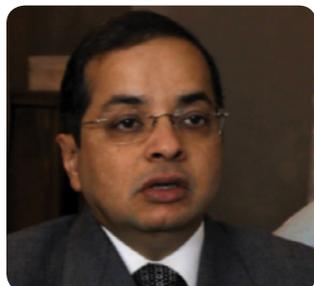

"Patients want two things... one is obvious relief from pain, but more importantly they hope that their replacement knee will last."

Dr Hemant Wakankar^{6*}
Orthopedic Surgical Consultant
Devchhaya Clinic
Maharashtra, India

Pre-clinical testing has demonstrated the potential for a femoral component made of PEEK-OPTIMA to offer equivalent performance to metal, while maintaining the benefits of a non-metallic solution.⁶ PEEK-OPTIMA exhibits properties closer to cortical bone than metal or ceramic materials, including its elastic modulus and density, and may offer a more natural solution to knee replacement.⁷⁻⁸ Used in approximately nine million implanted devices worldwide, including spinal cages and bone anchors, PEEK-OPTIMA has a strong clinical history with patients benefitting from this less stiff material.⁹



"Invibio is best placed to lead the development of a new innovative knee replacement solution because they have been successful in other areas and other arenas."

Dr Asit Shah^{11*}
Orthopedic Surgical Consultant
Englewood Orthopedic Associates
Englewood, NJ



"Technically it's very challenging but it's a very innovative concept and that's going to stimulate the team to rise to that challenge."

Professor John Fisher CBE^{*}
Director Institute of Medical & Biological Engineering
University of Leeds



"It will give designers a real opportunity in the future to change the way that we think about knee replacements."

Dr Jon Conroy^{12*}
Consultant Orthopedic Surgeon
Harrogate District General Hospital
North Yorkshire, England

Based on this pre-clinical work, a pathway for clinical trials to demonstrate PEEK-OPTIMA Knee's safety and clinical efficacy has been provided. These trials will assess post-operative knee performance and patient satisfaction levels. We are confident clinical trial results will show PEEK-OPTIMA Polymer Knee improves TKR outcomes and patient quality of life, while raising the bar for patient care. ▲

ABOUT THE AUTHOR

Ian Revie, PhD

Dr Ian Revie is Invibio Biomaterial Solutions' EU Developing Markets Marketing Leader responsible for the PEEK-OPTIMA Knee program. Prior to Invibio, Dr Revie was CEO for Activ4Life Healthcare Technologies where he patented a novel patient monitoring application focused on orthopedics. Previously, Dr Revie held progressively responsible roles in Johnson & Johnson's orthopedics business where he was credited for turning technical visions into commercial success. In 1995, Dr Revie received his PhD in biomechanical engineering from Queen's University, Belfast, Northern Ireland.



REFERENCES

1. Organization for Economic and Co-operation Development (OECD). (2015) *OECD Health Statistics*. Retrieved August 15, 2017, from <http://dx.doi.org/10.1787/health-data-en>.
2. Holy Cross Orthopedic Institute, The Leone Center for Orthopedic Care. (n.d.) *The Unhappy Total Knee Replacement Patient: Figuring out What Is Wrong*. Retrieved January 26, 2016, from <http://holycrossleonecenter.com/blog/unhappy-total-knee-replacement-patient-figuring-wrong/#sthash.FWSavHSE.dpuf>.
3. Bourne RB, et al. (January 2010) Patient Satisfaction after Total Knee Arthroplasty: Who is Satisfied and Who is Not? *Clinical Orthopaedics and Related Research*, 468:1, 57-63.
4. Scott CEH, Bugler KE, Clement ND, MacDonald D, Howie CR, and Biant LC. (June 25, 2012) Patient expectations of arthroplasty of the hip and knee. *The Bone and Joint Journal*. doi: 10.1302/0301-620X.94B7.28219.
5. Siebert W¹, Mai S, Kober R, Heeckt PF. (September 2002) Technique and first clinical results of robot-assisted total knee replacement. *The Knee*, 9:3, 173-80.
6. Since 2012, Dr. Hemant Wakankar has provided ad hoc consultancy services to Invibio Ltd.
7. Cowie RM, Briscoe A, Fisher J, and Jennings LM. (2016) *The Influence of Environmental Conditions on the Wear of an All-Polymer Knee Implant*. Poster #1913 presented at the 2016 Orthopaedic Research Society Annual Meeting, Orlando, FL, USA.
8. Choi K¹, Kuhn JL, Ciarelli MJ, and Goldstein SA. (1990) The elastic moduli of human subchondral, trabecular, and cortical bone tissue and the size-dependency of cortical bone modulus. *Journal of Biomechanics*, 23:11, 1103-1113.
9. Measured density of PEEK-OPTIMATM Natural lies within the accepted density of trabecular bone (1.02-1.20 according to Stendig-Lindberg et al. 1993).
10. Data on file with Invibio Biomaterial Solutions, Invibio Ltd.
11. Since 2012, Dr. Asit Shah, has provided ad hoc consultancy services to Invibio Ltd.
12. Since 2014, Dr. Jon Conroy, has provided ad hoc consultancy services to Invibio Ltd.

*The testimonials presented have been provided by practicing orthopedic surgeon(s) and researcher(s). Their view and experiences are their own and do not necessarily reflect those of others. "Invibio" disclaims any liabilities or loss in connection with the information herein.

Copyright ©2017 Invibio Ltd. INVIBIO™, PEEK-OPTIMA™, INVIBIO BIOMATERIAL SOLUTIONS™ are trademarks of Victrex plc or its group companies. All rights reserved.

**Invibio Ltd.**

Victrex Technology Centre
Hillhouse International
Thornton-Cleveleys
Lancashire
FY5 4QD, UK

Tel: +44 (0) 1253 898 000

FAX: +44 (0) 1253 898 001

Invibio Inc.

300 Conshohocken State Road
West Conshohocken, PA
19428
USA

Toll Free: 866-INVIBIO (468-4246)

Tel: (484) 342-6004

Fax: (484) 342-6005

**For further information please email us at info@invibio.com
or visit our website at:**

► Invibio.com

Victrex plc and/or its group companies ("Victrex plc") believes that the information in this document is an accurate description of the typical characteristics and/or uses of the product or products, but it is the customer's responsibility to thoroughly test the product in each specific application to determine its performance, efficacy, and safety for each end-use product, device or other application. Suggestions of uses should not be taken as inducements to infringe any particular patent. The information and data contained herein are based on information we believe reliable. Mention of a product in this document is not a guarantee of availability.

Victrex plc reserves the right to modify products, specifications and/or packaging as part of a continuous program of product development. Victrex plc makes no warranties, express or implied, including, without limitation, a warranty of fitness for a particular purpose or of intellectual property non-infringement, including, but not limited to patent non-infringement, which are expressly disclaimed, whether express or implied, in fact or by law.

Further, Victrex plc makes no warranty to your customers or agents, and has not authorized anyone to make any representation or warranty other than as provided above. Victrex plc shall in no event be liable for any general, indirect, special, consequential, punitive, incidental or similar damages, including without limitation, damages for harm to business, lost profits or lost savings, even if Victrex has been advised of the possibility of such damages regardless of the form of action.

Supporting information is available on request for all claims referenced in this document.

Copyright ©2017 Invibio Ltd. INVIBIO™, JUVORA™ PEEK-OPTIMA™, INVIBIO BIOMATERIAL SOLUTIONS™ are trademarks of Victrex plc or its group companies. All rights reserved.